


**WE RECOMMEND  
READING  
THE INSTRUCTION  
before use!**



the latest version  
Of the instructions  
you can always find on our Homepage:  
**[www.epads.at](http://www.epads.at)**

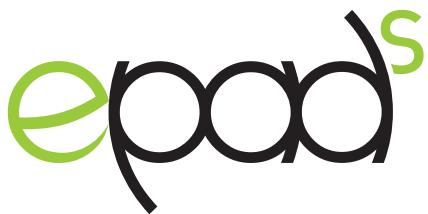

epad's

# **CONGRATULATION**

We are very happy that you have decided to purchase a the ePAD Sport and wish you a ing and wonderful feeling of wellness in using our product.

In ancient times the people used the positive effects of the magnetic field. The magnetic field is one of the oldest powers on our planet and has always been influential on the animal- and flora world. Also for us the earth magnetic field is vital. It gives us energy and protects us for cosmic radiation. The body regains the energy he loses through environmental influences in everyday life.

You have the possibility to use the ePAD Sport as your "private" magnetic field wherever you are.

## **ENERGY – ING- MOBILITY**

No matter where you are, after sport activities, on longer journeys, in the office, at home or on travel. The ePAD Sport should always be your constant companion, whenever you have the urge to .

The following pages of this instruction inform you about all the details and the easy use of the ePAD Sport.

We wish lots of fun and success with our product.

Your ePAD-Team

We work together with  
the Energy of the earth !

# 1. SIGNS AND SYMBOL



Read manual carefully



do not  
through into open fire



do not through into litter



rechargeable




do not throw into water



## SCOPE OF DELIVERY:


### ePAD Sport

- a) Controller
- b) Charger
- c) Charger cable
- d) ePAD
- e) Instruction

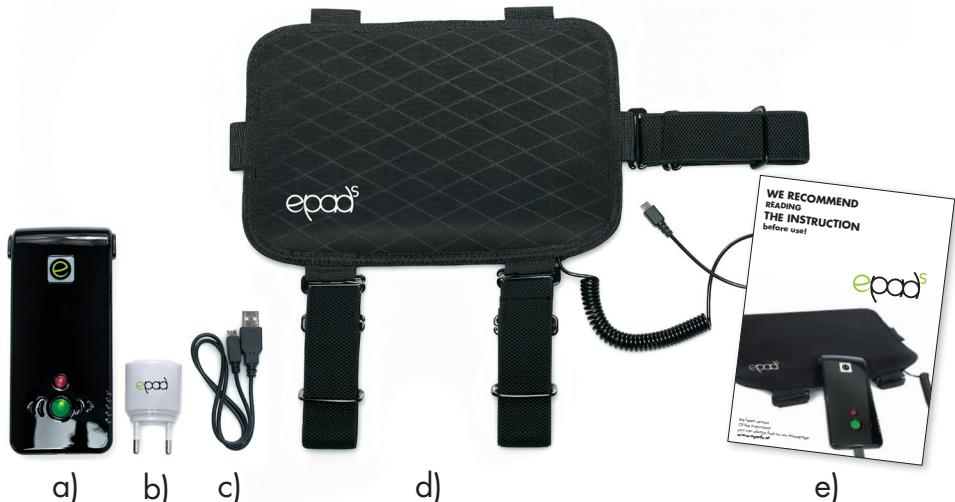

## 2. EXPLANATION OF THE DEVICE

The ePAD Sport is a portable and pulsating magnetic field Stimulation device, with which you can treat particular areas optimally or just boost up your energy.

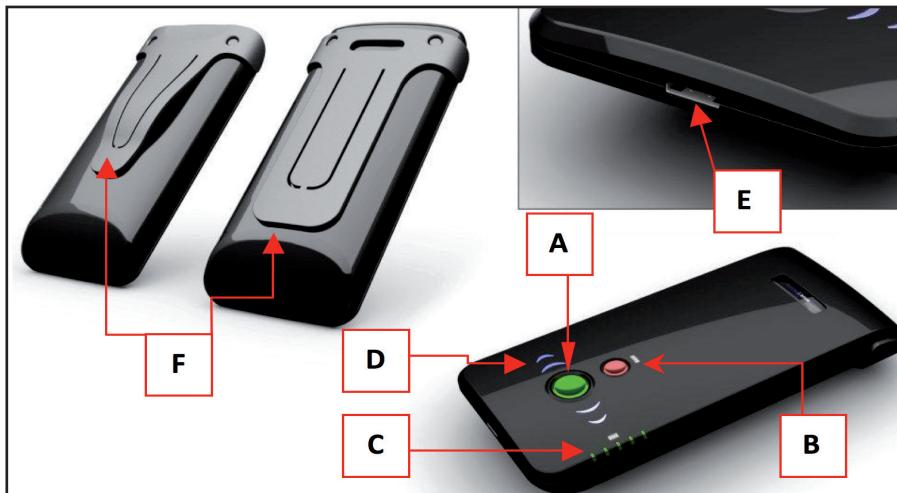

- A** on/ off button + program selection
- B** battery check (refer to page 8)
- C** battery charging display
- D** program display
- E** connection for ePAD Sport/connection of charging cable  
(NEW with USB connection cable)
- F** restraint clip

The controller must not get in contact with water.

## CONTRAINDICATIONS

**Persons** with electronic implants, pregnancy and persons with epileptic fits should not use the magnetic field

### **3. TECHNICAL DATA**

<b>IMPULSIVE FORM:</b>	Rectangular frequency
<b>FREQUENCY:</b>	3 hertz, 8 hertz, 1500 hertz
<b>POWER:</b>	<b>Rechargeable</b> Lithium-ion battery (4000mA, 3,6 volt), micro-USB To charge up the battery please use the charger enclosed.
<b>BATTERY POWER*:</b>	Depends on the frequency and temperature (at 20°C /68°F) approx 5-6 hours
<b>STORAGE TEMPERATURE:</b>	-20° C to +40° C / -4° F to 113° F
<b>INTERNAL MATERIAL:</b>	EVA-foam
<b>EXTERNAL MATERIAL:</b>	polyamide
<b>CONTROLLER:</b>	plastic

**IMPORTANT: recharge the controller continuous of a minimum of 8 hours for the first time use.**

**\*The battery power varies according to the surrounding temperature and program choice!**

### **4. WARRANTY**

MAS medical product trading company gives 2 (two) years of warranty on the ePAD Sport from the date of purchase to the first customer/ first buyer. It is only valid with prove of purchase. Please keep the receipt on a save place.

This warranty expands to cover the controller and also the applicator and covers all the defects, which are due to bad quality of the material and/or faulty installations.

## **NOT INCLUDED IN THE WARRANTY ARE:**

- Battery
- Damage caused by improper handling, accident or disregard to important hints
- Cracked or broken casing
- Yanked cable and/or damages on the plug
- Attempt of self-repair and/or own battery replacement

This guarantee does not affect the statutory rights of the consumers under the applicable national laws in the consumer's rights against the dealer arising from their sales/purchase contract.

## **NOTE!**

Service and repairs or a necessary exchange of battery may only be done by the manufacturer or by certified specialists. By non-compliance the warranty becomes invalid!

## **RETURN IN THE EVENT OF WARRANTY:**

In the event of warranty the device has to be sent back to MAS. If the device is returned under warranty, the buyer bears the cost of shipping to MAS, while MAS will take over the dispatch cost when returning the device back to the buyer.

## 5. CARE AND MAINTENANCE

The following instructions will help you to fulfil the warranty conditions.

- Clean the surface of the ePAD Sport regularly with a damp (not wet) cloth and let it air-dry.
- Do not put the ePAD Sport on source of heat (a radiator, heater or tiled stove etc.)
- Do not use any detergent or abrasive cleaner (for example chemical cleaning agents or steel wool).
- Never put the ePAD Sport in the washing machine or dryer
- Keep the ePAD Sport always on a cool and dry place
- Take out the USB stick carefully from the controller
- The controller must not get in contact with water.

## 6. DISPOSAL

To protect the environment the device including the battery should be disposed of as hazardous waste. Thank you!

## 7. INSTRUCTION FOR USE

### CONTROLLER

Take the controller and connect it to the cable of the applicator. Now press the red button (**B**) to test the battery power.

It shows you the current battery charged condition. The higher the amount of green glowing diodes (**C**) the better is the battery charged condition.



## HOW TO PUT ON THE RESTRAINT CLIP

Attach the elastic tension belt in the appropriate clips and adjust the length of the belts depending on the desired application.



## STARTING THE EPAD SPORT AND PROGRAM CHOICE:

Press the green button (**A**) **once (1x)**: now both waves (D) begin to flash **RED** every 4 seconds. Now you are in **program 1 (1500 hertz)**. This frequency corresponds to the acute inflammatory and energetic support in case of injuries.

Press the green button (**A**) **twice (2x)**: now both waves (D) flash **BLUE** every 4 seconds. Now you are in **program 2 (8 hertz)**. This represents the frequency of the organism in the sleeping phase and es the body.

Press the green button (A) three times (3x): now both waves (D) flash **GREEN** every 4 seconds. You are now in **program 3 (3 hertz)**. This represents the frequency of the organism in deep sleep; the body es completely and cramps disappear.

Press the green button (A) four times (4x): the device is turned off.

The device runs when the two waves (D) flash every 4 seconds at the same time. Now you can fix the controller on the belt clip (F). The program runs automatically for 25 minutes. Afterwards the device turns itself off automatically and the program is finished. You can stop a chosen program any time and by pressing the blue button (A) and change in another program or restart a program.

## TO ATTACH THE CONTROLLER PROPERLY



## EXAMPLES OF APPLICATION:



## TROUBLESHOOTING

### Error

No light shows

### cause

battery empty

### solution

connect the controller to the charger

### Error

Right wave doesn't light up

### cause

contact interruption

check the USB connection



epod's

## **Copyright 2011 - MAS medical trading company**

Instruction ePAD Sport; first edition; is valid from 1st of September 2011. All rights reserved. All form of duplication in whole or in extracts, or translation according to the Copyright law is not permitted unless approved by MAS in writing.

### **MAS medical trading company**

Hauptstrasse 50, A-8431 Gralla, [www.mas.at](http://www.mas.at)

# **Warranty card**

Please fill all the points in BLOCK LETTERS and attach the warranty card and prove of purchase when returning the product.

Serial number: .....

Date of purchase: .....

Name: .....

Address: .....

Short description of

fault: .....

Dealer's stamp: